

### Warnings

- ▶ For single use only. To avoid serious risk of infection, do not re-use or share the pen needle.
  - ▶ Keep it out of the children's reach. Not for use by children without adult supervision.
  - ▶ Consult your healthcare professional for assistance in choosing the right needle length, injection technique and site. Read the needle gauge and length on the seal to make sure you are using the right pen needle for the injection.
  - ▶ The pen needle is not recommended for people with visual impairment without the assistance of a sighted and trained individual.
  - ▶ Pen needle is sterile if the cap is unopened and undamaged. Do not use the pen needle if the individual packaging is damaged or if the seal is missing or if the needle appears to be bent or damaged.
  - ▶ To avoid accidental needle sticks, never handle someone else's pen needle.
  - ▶ Pen needles should be disposed safely in accordance with local regulations.
  - ▶ Follow the pen's Instructions For Use. Refer to the Instructions For Use from the insulin / pen injector manufacturer for specific storage and handling conditions. Seek immediate medical attention if you have any health concerns.
  - ▶ If you experience any serious incident that occurred in relation to the use of this product, please report it to the manufacturer and the competent authority of medical devices in your country.
- ▶ A serious incident means any incident that directly or indirectly led, might have led, or might lead to any of the following:
- (a) the death of a patient, user, or other people,
  - (b) the temporary or permanent serious deterioration of a patient's, user's or other person's state of health,
  - (c) a serious public health threat.

### Intended Use


FORAFINE disposable insulin pen needle is for single patient use only. It is intended for the subcutaneous injection of insulin with insulin pens.

### Operating Instructions

#### IMPORTANT





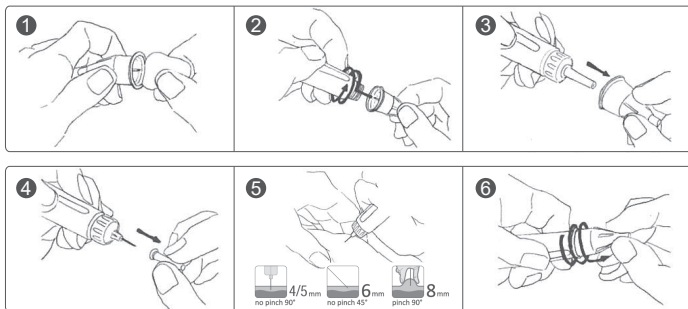







The injection site and technique should be recommended to you by your healthcare professional.


Choose and then disinfect the injection site. Change the injection site often to avoid tissue hardening. You can inject into the fatty tissue of the abdomen, buttocks, or the outer side or front of the upper thigh. Injections into the upper arm or the inner side of the thigh are not recommended, as the risk of injecting insulin into a muscle is high.




1. Remove the protective seal.
2. Press the needle straight on the pen and screw it clockwise until it is firmly attached.
3. Remove the outer protective cap.
4. Remove the inner protective cap.
5. Insert the needle on the injection site to deliver insulin.
6. After injection, replace the outer cap and unscrew the needle from the cap. Dispose of the needle safely.

## Symbol Information

| Symbol  | Referent                        | Symbol  | Referent  |
|---|---------------------------------|---|---|
|  | Consult instructions for use    |  | Manufacturer                                    |
|  | Caution                         |  | Do not use if package is damaged                |
|  | Use by / Expiry date            |  | Do not re-use                                   |
|   | Sterilized using ethylene oxide |  | Batch code                                      |
|   | CE mark                         |  | Model number                                    |
|   | Medical device                  |  | Authorised representative in the European Union |

 **ForaCare Suisse AG**  
Neugasse 55, CH9000, St. Gallen, Switzerland  
www.foracare.ch

 **MedNet EC-REP GmbH**  
Borkstraße 10, 48163 Münster,  
Germany

   ACS022    
  