

Easy Testing Steps:

Blood glucose test:



1 Insert test strip



2 Apply blood



3 Blood glucose data displayed

Blood pressure test:



1 Apply the cuff



2 Push the $\text{\textcircled{P}}$ button



3 Blood pressure and heart rate displayed

Specifications:

Device	FORA DUO ultima D40
Measurement Range	20 ~ 600 mg/dL (1.1~33.3 mmol/L)
Test Strip	FORA DUO ultima D40 test strips
Code	No coding required
User Selection	4-user option
Sample Volume	0.5 μ L
Reaction Time	5 seconds
Coefficient of Variation	Less than 5%
Pressure Range	0 ~ 300 mmHg
Heart Rate Range	40 ~ 199 beats/min
Memory Capacity	864 results with date & time
Dimension	147 x 105 x 80 mm
Weight	400 g (batteries excluded)
Power Source	AA*4 + 6V Adapter or Li-ion Battery
Data Output	USB or Bluetooth or GPRS connectivity

ForaCare Suisse AG
www.foracare.ch

Neugasse 55
9000 St. Gallen, Switzerland
Tel: +41-(0)71-220-1001 Fax: +41-(0)71-220-1075

FORA®

DUO ultima D40

2-in-1 Blood Glucose plus Blood Pressure Monitoring System

Best for Dia-hyper-besity



Automated Averaging Function




Talking-Voice Guided



PC Link

Why 2-in-1? Diabetes & Hypertension

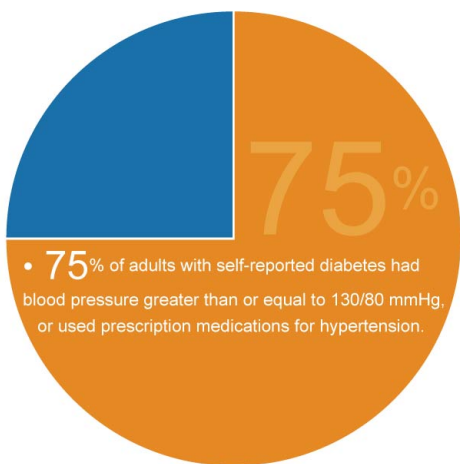

Compared to the general population, people with diabetes are confronted by a twice to fourfold increased risk of suffering from cardiovascular diseases.

A diabetes patient with hypertension:

- triples the risk of coronary artery disease
- doubles the overall mortality rate and risk of stroke

Furthermore, hypertension also accelerates the deterioration of diabetic complications including diabetic nephropathy, retinopathy and neuropathy. Therefore, self monitoring of these two crucial signs is of urgent necessity. The 2-in-1 monitoring system by FORA has made measurements of BG and BP accurate, safe and easy, all within ONE device.

Healthy lifestyle, regular measurement and prudent management of blood glucose plus blood pressure values are equally crucial for people with or without diabetes and hypertension. Please consult medical professionals for your specific blood glucose and blood pressure safe zones, and start an accurate and easy self monitoring with FORA today!

Reference:
Diabetes, hypertension, and cardiovascular disease: an update. Hypertension 37:1053-1059, 2001
Diabetes and atherosclerosis: epidemiology, pathophysiology, and management. JAMA 287:2570-2581, 2002

DUO ultima D40

2-in-1 Blood Glucose plus Blood Pressure Monitoring System

AVG 3 times average

Automated Averaging Function

2in1 Blood Glucose plus Blood Pressure

ONE meter / TWO tests

Strip Ejection

Avoid infection

4-user selection

4-user Selection

Comfortable Easy Use

Auto inflation cuff with Smart Pumping


PC Link

Data transmission for better management




Less Pain with **AST**

Alternative Site Testing
Blood sample can be obtained from palms, forearms, upper arms, calves, and thighs -- with proven equivalent accuracy to fingertip testing, thus reducing pain on fingertips.

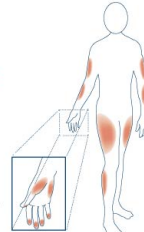

Data Management:

Testing data management is important for controlling both blood glucose and blood pressure. FORA Telehealth System assists users to transfer, organize and analyze data efficiently and wisely.

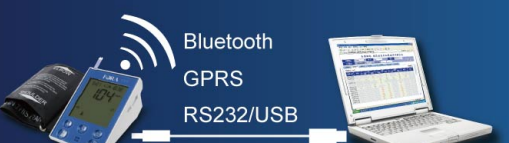

Please learn more from www.foracare.com

